


Upright Ergometer 8.5UE



The 8.5UE upper body ergometer offers a very low starting resistance to accommodate various ranges of motion capabilities, ensuring gentle and low-impact movements for the upper limbs, arms, and shoulders. Additionally, the adjustable crank length can be tailored to meet each patient's specific range of motion requirements. The 8.5UE is thoughtfully designed with the user in mind. It is feature-rich and ADA-compliant. With a fully removable seat, it facilitates seamless direct wheelchair access. Users can choose from three distinct workload options, including Constant Power, Isokinetic, and Constant Resistance, making it adaptable to a wide range of rehabilitation needs. Moreover, the system boasts a suite of neurological programs and features tailored to rehabilitation. The 10-inch touchscreen display is equipped with multi-connectivity protocols and provides access to cloud analytics, making the 8.5UE an ideal tool for both therapists and patients, promoting enhanced therapy and rehabilitation outcomes.

- ADA-compliant with direct wheelchair access
- Adjustable pedal crank to fit any range of motion
- 50 levels of forward and reverse pedaling resistance
- Large round rotating hand grips
- Choose from a variety of preset programs or create a custom program
- Supportive high-density foam seat with fore, aft, and recline angle adjustments
- Accommodates users between 4'10" - 6'6" and up to 440 lb

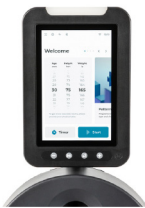

Intuitive, easy-to-read display



Adjustable crank arm length for patients with limited mobility



Multi-function handle with rotating hand grips, remote shortcut keys and arm height adjustment lever.

Specifications

| | |
|-------------------------|--|
| Console | 10.1" Tablet Console, USB for data transfer software & cardiac stress test, 360 degree horizontal adjustment |
| Resistance System | Electromagnetic resistance system, 50 levels |
| Heart Rate Receiver | Bluetooth, ANT+ |
| Programs | Manual / Timer / Hill / Plateau / HIIT / Facility / Heart Rate Target / Calories Target / METs Target / Wingate Test / Steep Ramp Test / Conconi Test / Power Distribution |
| Power | 100-240V Cord/Outlet, 3V CR1220 Battery |
| Crank Arm | 12 Levels |
| Seat | High Density Foam, 360-degree Swivel, Fore/Aft/Seat Back Angle Adjustments, Removable for Wheelchair Access |
| Workload Range | 5 to 1000 Watts |
| Recommended User Height | 4'10" - 6'6" |
| Max User Weight | 440 lb |
| Dimensions | 67.2" x 29.9" x 65.9" |
| Handrail | Round rotating hand grips |
| Product Weight | 328.5 lb |
| Compliance | MDD 93/42/EEC Class Im cTUVus |
| Warranty | 5 Years parts, 1 Year labor |